




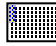




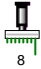

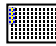



1	Set Variable	VolOrangeG "Volume to check pipetting precision?" "Volume to check pipetting precision?", 0.25 - 10
2	Set Variable	VolWater = 200-VolOrangeG
3	Wash Tips	 25 + 0.1 ml
4	Wash Tips	 4 + .1 ml
5	Wash Tips	 .2 + .1 ml
6	Begin Loop	3 times "Loop1"
7	Aspirate	 VolWater*4 µl Water free dispense "Labware9" (Col. 1, Rows 3-6)
8	Begin Loop	4 times "Loop"
9	Dispense	  VolWater µl Water free dispense "Labware18" (Col. 1, Rows 1-4) , 2 options
10	End Loop	"Loop"
11	Wash Tips	 5 + 4.0 ml
12	End Loop	"Loop1"
13	Begin Loop	12 times "Loop3"
14	Aspirate	 VolOrangeG µl D-BSSE Te-PS Wet Contact "Labware5" (Col. 1, Row 1)
15	Aspirate	 VolOrangeG µl D-BSSE Te-PS Wet Contact "Labware5" (Col. 1, Row 1)
16	Aspirate	 VolOrangeG µl D-BSSE Te-PS Wet Contact "Labware5" (Col. 1, Row 1)
17	Aspirate	 VolOrangeG µl D-BSSE Te-PS Wet Contact "Labware5" (Col. 1, Row 1)
18	Dispense	  VolOrangeG µl D-BSSE Te-PS Wet Contact "Labware18" (Col. 1, Rows 1-4) , 1 option
19	Wash Tips	 5 + 4.0 ml
20	End Loop	"Loop3"
21	Transfer Labware	Source: Grid '10,' Site '1'; Destination: Grid '10', Site '2'; Narrow (ROMA 2); Move to home position
22	H+P Shaker	HP ShakeForTime(*271 8*30*30*30*30 2 *30 1 *30 1,5*30*30 255*27)
23	Start Timer	1

24	Wait for Timer	Timer 1 : 120 sec
25	H+P Shaker	HP ShakeForTime(*271 8*30*30*30*30 2 *30 1 *30 1,5*30*30 255*27)
26	Start Timer	1
27	Wait for Timer	Timer 1 : 120 sec
28	H+P Shaker	HP ShakeForTime(*271 8*30*30*30*30 2 *30 1 *30 1,5*30*30 255*27)
29	InfiniteM200	Open()
30	Transfer Labware	Source: Grid '10,' Site '2'; Destination: Grid '61', Site '1'; Narrow (ROMA 2); Move to home position
31	Transfer Labware	Source: Grid '61,' Site '1'; Destination: Grid '68', Site '1'; Wide (ROMA 2); Move to home position
32	InfiniteM200	Measure(C:\Program Files\TECAN\EVOware \output\Infinite\<Y YYYY-MM-DD HH-mm-ss> .xls <TecanFile xml ns:xsiequal; quot e; http://www.w3.or g/2001/XMLSchema-in stancequote; xsi: schemaLocation_equa l; quote; tecan.at .schema.documents M ain.xsdquote; fil eformatequal; quo te; Tecan.At.Measur ementquote; filev ersionequal; quot e; 2.0quote; xmln sequal; quote; te can.at.schema.docum entsquote; ><FileI nfo typeequal; qu ote; quote; instr umentequal; quote ; infinite 200Proq uote; version_equa l; quote; quote; createdFromequal; quote; installqu ote; createdAtequa l; quote; 2011-03 -04T09:52:03.9375Z quote; createdWith equal; quote; Tec an.At.XFluor.Reader Editor.XFluorReader Editorquote; desc riptionequal; quo te; quote; /><Tec anMeasurement idequal; quote; 1quot e; classequ al; q uote; quote; ><Mea surementManualCycle idequal; quote; 2quote; numbereq ual; quote; 1quot e; typeequal; qu ote; Standardquote ; ><CyclePlate ide qual; quote; 3quo te; fileequal; g uote; GRE96ftquote ; plateWithCover equal; quote; False quote; ><PlateRang e idequal; quote; 0quote; rangeeq ual; quote; A1til de; H12quote; aut oequal; quote; tr uequote; ><Measure mentAbsorbance ide qual; quote; 4quo te; modeequal; q uote; Normalquote; typeequal; quot e; quote; nameeq ual; quote; ABSqu ote; longnameequa l; quote; quote; descriptionequal; quote; quote; >< Well idequal; quo te; 5quote; auto equal; quote; true quote; ><Measureme ntReading idequal; quote; 6quote; nameequal; quote; quote; beamDiamet erequal; quote; 700quote; beamGri dTypeequal; quote ; Singlequote; be amGridSizeequal; quote; 1quote; be amEdgeDistanceequa l; quote; autoquo te; ><ReadingLabel idequal; quote; 7 quote; nameequal ; quote; Label1qu ote; scanTypeequa l; quote; ScanFixe dquote; refIDequa l; quote; 0quote ; ><ReadingSettings numberequal; quo te; 10quote; rate equal; quote; 250 00quote; /><Readi ngTime integrationT imeequal; quote; 0quote; lagTimee qual; quote; 0quote; readDelayequa l; quote; 100000q uote; flashequal; quote; 0quote; darkequal; quote; 0quote; /><Readi ngFilter idequal; quote; 0quote; t ypeequal; quote; Exquote; waveleng thequal; quote; 4 920quote; bandwid thequal; quote; 9 0quote; attenuati onequal; quote; 0 quote; usageequa l; quote; ABSquote ; /></ReadingLabe l></MeasurementRead ing></Well><Well id equal; quote; 8q uote; autoequal; quote; truequote; ><MeasurementReadi ng idequal; quote; 9quote; nameeq ual; quote; quote ; beamDiameter_equa l; quote; 700quo te; beamGridTypee qual; quote; Singl equote; beamGridS izeequal; quote; 1quote; beamEdgeD istanceequal; quo te; autoquote; ><R eadingLabel idequa l; quote; 10quote ; nameequal; quo te; Ref1_Label1quo te; scanTypeequal ; quote; ScanFixed quote; refIDequa l; quote; 1quote; ><ReadingSettings numberequal; quot e; 10quote; rate equal; quote; 2500 0quote; /><Readin gTime integrationTi meequal; quote; 0 quote; lagTimeeq ual; quote; 0quot e; readDelayequal ; quote; 100000qu ote; flashequal; quote; 0quote; d arkequal; quote; 0quote; /><Readin gFilter idequal; quote; 0quote; ty peequal; quote; E xquote; wavelengt hequal; quote; 62 00quote; bandwidt hequal; quote; 90 quote; attenuatio nequal; quote; 0 quote; usageequal ; quote; ABSquote ; /></ReadingLabel ></MeasurementReadi ng></Well><CustomDa ta idequal; quote ; 11quote; /></Me asurementAbsorbance ></PlateRange></Cyc lePlate></Measureme ntManualCycle><Meas urementInfo idequa l; quote; 0quote; descriptionequal ; quote; quote; > <ScriptTemplateSett ings idequal; quo te; 0quote; ><Scri ptTemplateGeneralSe ttings idequal; q uote; 0quote; Tit leequal; quote; _quote; Groupequal ; quote; quote; Infoequal; quote; quote; Imageequa l; quote; quote; /><ScriptTemplate DescriptionSettings idequal; quote; 0quote; Internal equal; quote; quo te; Externalequal ; quote; quote; IsExternalequal; quote; Falsequote; /></ScriptTemplat eSettings></Measure mentInfo></TecanMea surement></TecanFil e>)
33	Transfer Labware	Source: Grid '68,' Site '1'; Destination: Grid '61', Site '1'; Wide (ROMA 2); Move to home position
34	InfiniteM200	Close()

Transfer Labware

Source: Grid '61,' Site '1'; Destination: Grid '10', Site '1'; Narrow (ROMA 2); Move to home position

